

AERYS®



AERYS® bed with ABÉLIA II boards, Light Beech finish, RAL1001 handrails.




THANKS TO ITS UNDERSTATED STYLE, A REFRESHING CHANGE FROM THE USUAL MEDICAL ENVIRONMENT, AND ITS HOST OF FUNCTIONS, THE AERYS BED IS SUITABLE FOR A RANGE OF DIFFERENT PATIENT CONDITIONS AND IMPROVES ERGONOMICS FOR CAREGIVERS.

THE AERYS BED HAS A NEW EMBEDDED ELECTRONIC MANAGEMENT SOFTWARE PROGRAM IN ITS CONTROL BOX. THE AERYS BED COMES STANDARD WITH:

- ▶ Possibility of the variable height data to customize and facilitate getting out of bed for patients.
- ▶ One-touch access to the chair function.
- ▶ A stop in the -5° Trendelenburg position to reposition the patient.
- ▶ A Trendelenburg position at -12° for emergencies.

● If you choose the option (f) with the multi-directional head trolley, a beep and a 3-second stop will tell you that the position has switched from braked to un-braked.

You can also choose the epoxy coating color of your bed from our color chart.

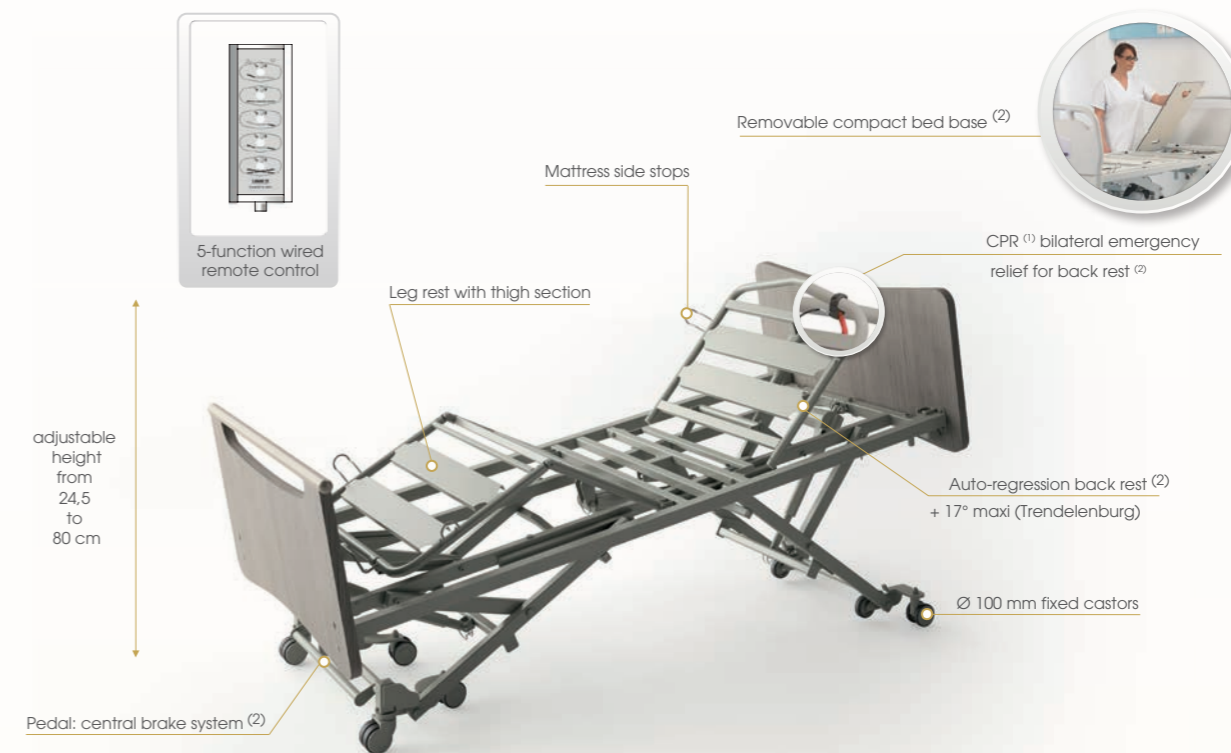

(1) CPR - Cardio Pulmonary Resuscitation

(2) OPTIONS AND ACCESSORIES

- ▶ DESIGN castors Ø 75mm
- ▶ Centralized braking on foot
- ▶ CPR electric back rest
- ▶ Auto-regression back rest
- ▶ Removable compact base board
- ▶ 4-castor multi-directional movement system Ø 75mm castors at the head
- ▶ EVOLUTION, ACTIV & MOBILITY AID SAM range
- ▶ Backup battery
- ▶ Remote control on flexible arm
- ▶ Caregiver's wired control console
- ▶ Infrared remote control
- ▶ Bed base extension
- ▶ Wall bumper
- ▶ Different epoxy coating colors to choose from

DESCRIPTION	DIMENSIONS	ADJUSTABLE HEIGHT	BR	LEG REST	LIFTING CAPACITY	CASTORS	STANDARDS	MOTORIZATION
AERYS								
3 electrical functions without knee break without rtr position 3 electrical functions with knee break rtr/tr position	90 x 200 cm Bed base non-divisible with welded metal slats	24.5 cm to 80 cm	70°	thigh section 26° and 16° feet adjustment by manual crank	SWL (1) 200 kg Patient weight up to 165 kg Accessories 15 kg Mattress 20 kg (Environment 3, 4 with no RTR/TR function and 5)	Head trolley with fixed metal castors Ø 100 mm Foot trolley with multi-directional castors and individual braking Ø 75 mm	IEC 60601 IEC 60601-2-52 compatible with environments 3	LINAK IP54

Handset, electric and mechanical parts: 5-year warranty*
Boards and accessories: 1-year warranty*
*Under normal conditions of use and excluding wearing parts - please see instructions of use

(1) Safe working load

Photos not contractually binding

FP000064/01.2018



WINNCARE Group - 4, Le Pas du Château - 85670 Saint-Paul-Mont-Penit - France
Tél : +33 (0)2 51 98 55 64 - Fax : +33 (0)2 51 98 59 07 - E-mail : info@medicatlantique.fr - www.winnicare.fr
In the interest of constantly improving its products, the WINNCARE Group reserves the rights to modify the technical specifications of any product without prior notice
Class I medical device in accordance with European Directive 93/42/EEC

